
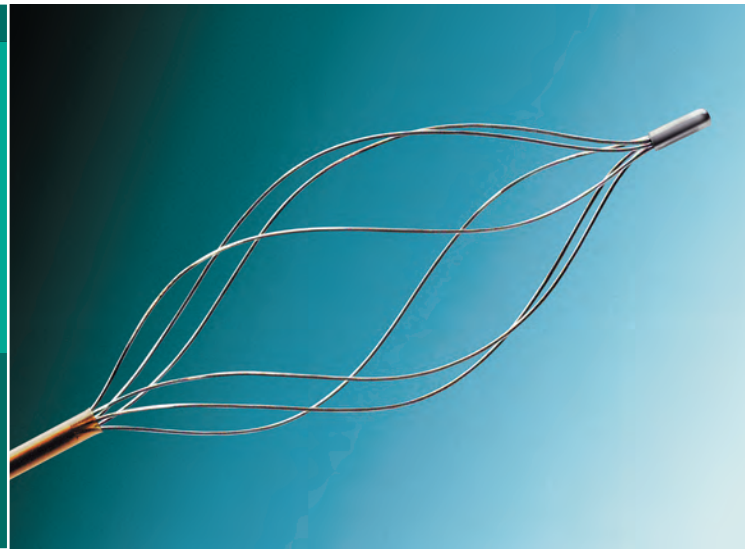


When the challenge
is efficient capture
of calculi ...

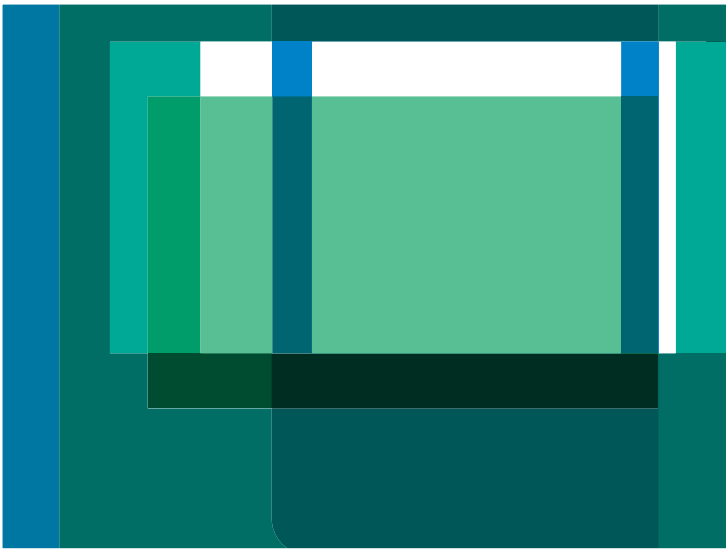

EXPAND²¹²®
HELICAL STONE BASKET

the possibilities
with an innovative
solution.

The latest innovation in stone basket design and performance

Take advantage of the exceptional performance characteristics of our innovative technology for calculi removal: the EXPAND²¹²® Helical Stone Basket. When the choice is for a helical wire design, the EXPAND²¹²® basket combines all the benefits essential to an efficient stone removal procedure: strength, durability and an unsurpassed capture rate.*




New Paired Wire Design




- Patented “2-1-2-1” design features alternating paired and single stainless steel wires for secure grasp of calculi and stone fragments.
- Kink-resistant, torque-stable Nitinol drive wire and proprietary sheath material promote effective navigation through tortuous anatomy.
- Round wire construction promotes atraumatic manipulation of calculi because there are no sharp edges.

Exceptional Expansion Force

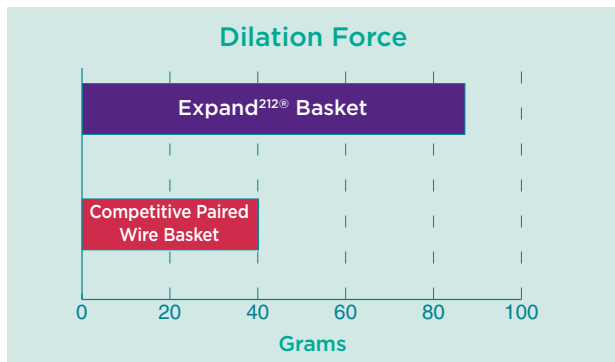

- Laboratory tests showed that the EXPAND²¹²® basket exhibited more than twice the dilation force of the leading competitive paired wire basket.*
- Exceptional dilation force optimizes performance by helping the EXPAND²¹²® basket retain its shape even after multiple passes.

Remarkable Capture Ability

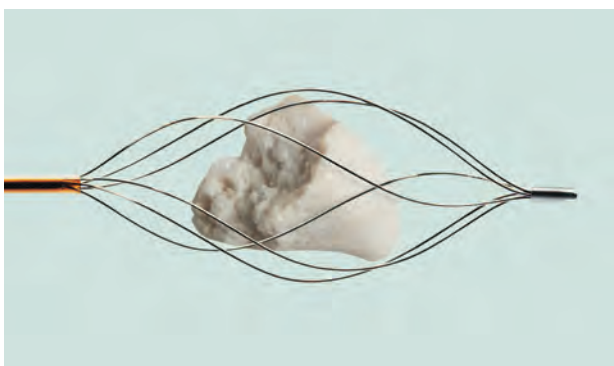

- Exclusive wire design minimizes lateral movement, permitting more definitive stone manipulation.
- EXPAND²¹²® basket opens effectively in the undilated ureter, allowing for the efficient capture of calculi and small stone fragments.

Excellent Control



- The EXPAND²¹²® basket features the exclusive Bard® Platinum Class™ handle for ergonomic comfort, tactile “feel” and control.
- The Nitinol drive wire provides torque stability for 1:1 rotation.
- Luer clip permits simple removal of handle from basket.

EXPAND²¹²®

HELICAL STONE BASKET

the possibilities with an innovative solution.

Innovative Design

- Patented paired wire design in “2-1-2-1” configuration helps assure effective capture of calculi and stone fragments.
- Nitinol drive wire and proprietary sheath material promote smooth, atraumatic insertion and navigation through delicate tissues.

Exceptional Performance

- Laboratory tests demonstrated an exceptional expansion force when compared with the leading competitive paired wire basket.*
- Bard® Platinum Class™ handle delivers ergonomic comfort, tactile feel and excellent control throughout the procedure.

Reorder Number	Basket Diameter (mm)	Length (cm)	Shaft Size (Fr.)	Filiform Tip
040490	11	90	3	No
040415	11	115	3	No
041490	11	90	3	5cm
041415	11	115	3	5cm

For additional information contact your local Bard representative or call 1.888.367.2273.



Where Knowledge Leads™



C. R. Bard, Inc. Covington, GA 30014 888.367.2273 www.bardurological.com

Bard and Expand 2121 are registered trademarks of C. R. Bard, Inc. or an affiliate. Platinum Class and Where Knowledge Leads are trademarks of C. R. Bard, Inc. or an affiliate. © 2008, C. R. Bard, Inc. Printed in U.S.A. * Data on file, C. R. Bard, Inc.

Please consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions, and directions for use.

0803-28 R04/08 THP P4/08 5M

*Data on file, C. R. Bard, Inc.